







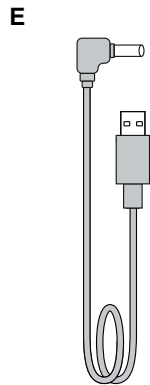
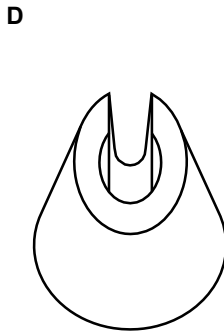
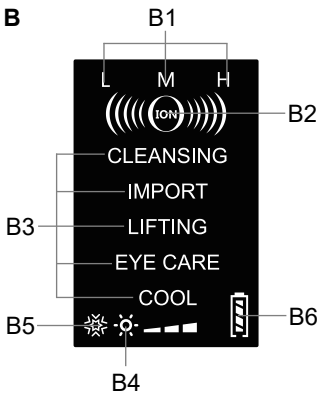
<i>Galvanická žehlička na pleť</i> • NÁVOD K OBSLUZE	CZ	5-10
<i>Galvanická žehlička na pleť</i> • NÁVOD NA OBSLUHU	SK	11-16
<i>Galvanic skin iron</i> • USER MANUAL	EN	17-22
<i>Galvanikus bőrvasaló</i> • KEZELÉSI ÚTMUTATÓ	HU	23-28
<i>Żelazko galwaniczne do cery</i> • INSTRUKCJA OBSŁUGI	PL	29-34
<i>Galvanisches Hautbügeleisen</i> • BEDIENUNGSANLEITUNG	DE	35-41

Fenité

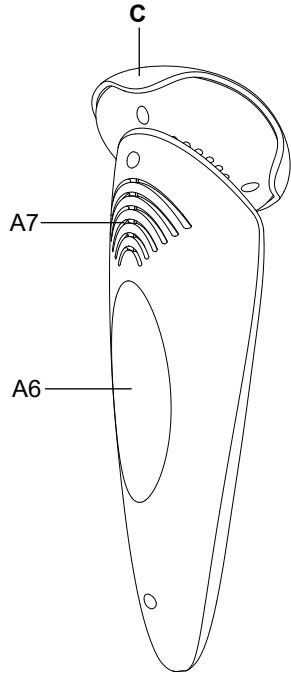1

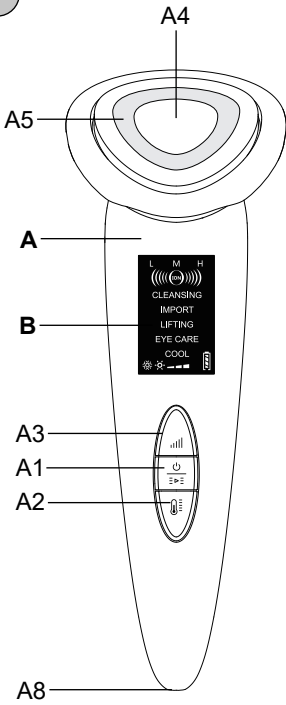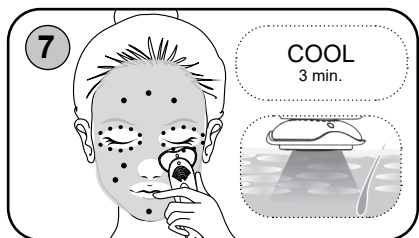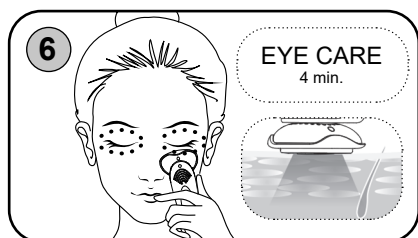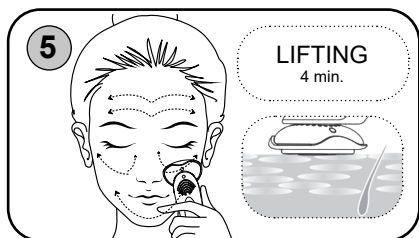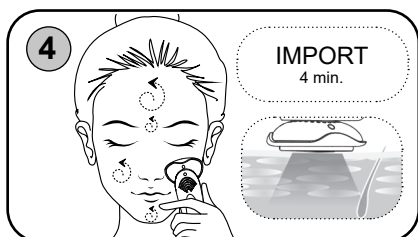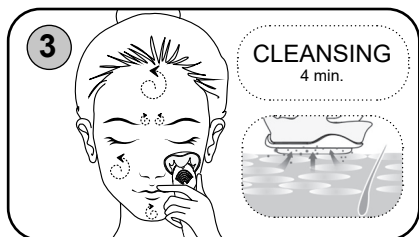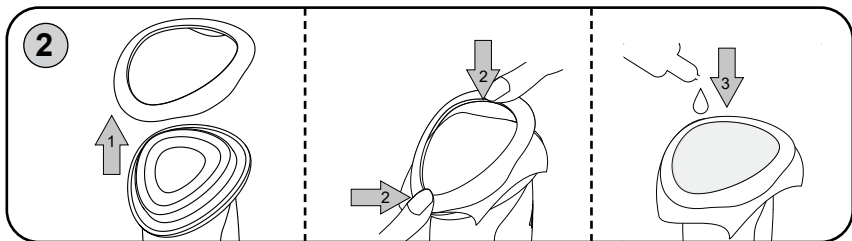CZ	I. BEZPEČNOSTNÍ UPOZORNĚNÍ	5
	II. POPIS SPOTŘEBIČE (obr. 1)	7
	III. POKYNY K OBSLUZE	7
	IV. ÚDRŽBA	9
	V. EKOLOGIE	10
	VI. TECHNICKÁ DATA	10
SK	I. BEZPEČNOSTNÉ UPOZORNENIA	11
	II. OPIS SPOTREBIČA (obr. 1)	13
	III. POKYNY K OBSLUZE	13
	IV. ÚDRŽBA	15
	V. EKOLÓGIA	16
	VI. TECHNICKÉ ÚDAJE	16
EN	I. SAFETY WARNING	17
	II. DESCRIPTION OF THE APPLIANCE (Fig. 1)	19
	III. INSTRUCTIONS FOR USE	19
	IV. MAINTENANCE	21
	V. ENVIRONMENT	22
	VI. TECHNICAL DATA	22
HU	I. BIZTONSÁGI RENDELKEZÉSEK	23
	II. A KÉSZÜLÉK LEÍRÁSA (1. ábra)	25
	III. HASZNÁLATI UTASÍTÁSOK	25
	IV. KARBANTARTÁS	28
	V. ÖKOLÓGIA	28
	VI. MŰSZAKI ADATOK	28
PL	I. OSTRZEŻENIA DOTYCZĄCE BEZPIECZEŃSTWA	29
	II. OPIS URZĄDZENIA (rys. 1)	31
	III. INSTRUKCJE DOT. OBSŁUGI	31
	IV. KONSERWACJA	34
	V. EKOLOGIA	34
	VI. DANE TECHNICZNE	34
DE	I. SICHERHEITSAUWEISUNGEN	35
	II. BESCHREIBUNG DES GERÄTES (Abb. 1)	37
	III. BEDIENUNGSANWEISUNGEN	38
	IV. WARTUNG	40
	V. UMWELTSCHUTZ	40
	VI. TECHNISCHE DATEN	40

eta 6353

Fenité

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INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please keep these instructions including the warranty, the receipt and, if possible, the carton box with the internal packing.

I. SAFETY WARNING

- Please consider the instructions herein as a part of the appliance and pass them over to the next user of the appliance.
- The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in beauty parlours or for any other commercial utilisation !
- This appliance may only be used by children aged 8 years or older, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never immerse the appliance in water!
- When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.
- Never use the appliance if it does not correctly operate, fell on the ground or suffered any damage.
- Protect the appliance from dust, chemicals, excessive humidity, intensive temperature changes, and do not expose it to direct sunlight.
- Before use, take off all metal items and jewellery (e.g., glasses, watch, earrings, bracelets, rings, etc.).
- Do not insert the appliance or its accessories in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- Only use the appliance with the accessories designed for this type. Use of unapproved accessories may cause injury.
- This appliance contains batteries that are not replaceable.
- This appliance has a heated surface, and the persons with impaired heat perception must be careful in using it.
- **WARNING:** In the case of incorrect use of the appliance that contradicts the operating instructions, there is a risk of injury.
- Texts in a foreign language and images provided on the packaging or on the product, if any, are translated and explained at the end of this language version.
- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) and is not responsible according to the appliance warranty in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.
- Do not use the appliance if your skin is oversensitive, damaged or exematic (e.g. in the case of fresh scars, small injuries, burns, erythema, swollen cheeks, skin diseases...).
- If you perceive pain/ unpleasant feelings or if any other health problems occur, switch off the appliance immediately and consult with your physician or a dermatologist about its use.
- Do not use this appliance if you have metallic implants, metallic dental braces, piercing or an implanted pacemaker.
- Consult the use of this appliance with your physician if you have any of the following contraindications: pregnancy, breastfeeding, allergy, blood coagulation disorder, chronic diseases or if you take sensitive medical preparations sensitive to light.

- The appliance is not intended to cure any skin disease symptoms. Please consult use of the appliance with your doctor, where necessary.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.

SPECIAL SAFETY INSTRUCTIONS

- Children must be kept away from the appliance and its adapter.
- To power the appliance, only apply safe small voltage in accordance with the type plate.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output, PC with USB output with output voltage of 5.0 V, and output current of at least 1000 mA). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.
- Never use the appliance if the adapter is damaged, if the appliance does not operate properly, fell on the ground and suffered any damage, or fell into water. In such case bring the appliance to the specialised service centre to verify the appliance safety and proper operation.
- For additional protection we recommend installing a protective switch (RCD) in the electrical circuitry of the bathroom, with the nominal actuating current not exceeding 30 mA. Ask an inspecting technician or an electrician for advice.
- The appliance and the connected adapter must not be used at places where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out. Then bring the appliance to a specialised service centre to check its safety and proper functionality.
- If the appliance is used in a bathroom, it must be disconnected from the mains after use by unplugging the power cord from the wall socket, because the appliance poses a risk if near to water, even when switched off.
- Always unplug the appliance from charger if you leave it unattended and cleaning.
- The mains adapter supplied with the appliance is only designated for this appliance.
- Never plug the adapter in the electrical socket or unplug it with wet hands.
- Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water!
- The appliance gets hot when charging, which is absolutely normal.
- Charge the battery at a temperature close to the room temperature.

- Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C, adversely affects the accumulator service life.
- Prevent hanging the power adapter cord loosely over desktop edge where children may reach it.
- The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.

Note

ETA a.s. is not a licensed provider of health and cosmetic care and does not have expertise in the field of treatment, diagnostics of the health condition and related effects or suitability for the use of appliances that may affect the health condition in any way. During any exercise, massage therapy or cosmetic intervention, there is a theoretical possibility of physical injury or, for example, worsening of the symptoms of existing diseases. By using this appliance, you agree that you use it at your own risk and assume all risks of injury or deterioration. You also agree to waive any claims (whether known or unknown) against ETA a.s. arising from improper use of this appliance or through negligence. Before using any appliances that may affect your health, first consult your doctor about the suitability of their use.

II. DESCRIPTION OF THE APPLIANCE (Fig. 1)

A – Appliance body

- A1 – ON/OFF button + operating mode setting button
- A2 – temperature setting button
- A3 – galvanisation intensity adjusting button
- A4 – application head
- A5 – LED (blue/red) for photon therapy
- A6 – contact area for correct galvanisation
- A7 – ventilation grille
- A8 – USB charging socket

B – LED display

- B1 – galvanisation intensity (L - low / M - medium / H - high)
- B2 – ionisation indicator
- B3 – operation modes (CLEANSING / IMPORT / LIFTING / EYE CARE / COOL)
- B4 – heating indicator (3 heating levels)
- B5 – cooling indicator
- B6 – Battery level indicator

C – Cloth fixing ring

D – Storage

E – USB cable *

* USB cable: A source with USB terminal may be connected to USB cable (for example AC power adapter, backup power source, etc.).

III. INSTRUCTIONS FOR USE

Before using the appliance, its battery must be charged first. Charging is possible with the enclosed USB cable from the mains, i.e. an AC adapter, or a backup power source.

CHARGING

- Before the first use or if the appliance was not used for a long time, charge it completely.
- The standard charging time is **approximately 2,5 hours**. and the fully charged accumulator allows for the operating time of **up to 5 hours**. The operating time decreases depending on the appliance loading and due to the accumulator ageing and wear and tear.

WARNING

Do not use the USB ports on other devices (e.g. computer, laptop, tablet, cell phone etc.) to charge the appliance. These devices do not have to be sized for the current consumption of the appliance and they might be damaged.

1. Connect the appliance to the power source (electric charger) using USB cable **E**. The appliance can also be used during the charging process.
2. Do not exceed the designated charging time!

Note

Battery level and charging state are displayed by indicator **B6**. The number of lighting segments increases during the charging process. When the charging process has been completed, charging indicator will light continuously.

Optimise the battery service life

- Do not leave USB cable connected to the appliance.
- During charging or regular operation, the appliance heats up, which is a normal phenomenon, and not a reason for a complaint.

CONTROL**Note**

In order to provide for proper galvanisation, area **A6** must be touched with your hand/fingers. Consequently, the appliance is intended exclusively for personal care.

On/off button + setting operation modes (A1)

- By long pushing of button **A1** (approx. for 2 sec.) you switch the appliance on/off. Upon switching on and off an acoustic signal sounds.
- Press short **A1** button to switch over the modes. – The operation mode you have set is shown on display **B** by lighting up its name (CLEANSING / IMPORT / LIFTING / EYE CARE / COOL).

Temperature setting button (A2)

- Briefly press button **A2** to adjust the temperature (3 heating levels are available). Set temperature is indicated by individual segments on display **B** (next to symbol **B4**).
- The temperature can be adjusted in operation modes CLEANSING, IMPORT a LIFTING.

Galvanisation intensity button (A3)

- Briefly press button **A3** to adjust galvanisation intensity (3 intensity levels). Set intensity level is visualised by illumination of the relevant symbol **B1** (L - low / M - medium / H - high) on display **B**.
- Galvanisation intensity can be adjusted in operation modes CLEANSING, IMPORT and LIFTING.
- In EYE CARE mode, the red light can be switched on/off by brief pressing button **A3**.

OPERATION MODES AND THEIR USE**CLEANSING**

This mode is intended for efficient and deep skin cleansing. It combines positive ionisation, heating and vibrations and allows removal of dead skin cells, sebum, impurities, traces of makeup and more.

Procedure (Fig. 3):

1. Place the cotton cloth on the application head **A4** and apply the makeup remover (Fig. 2).
2. Press **A1** to switch on the appliance, select the CLEANSING mode and take the appliance in your hand, touching area **A6** with your hand/fingers.
3. Where necessary and considering sensitivity of the skin, use button **A2** to adjust temperature and button **A3** to adjust galvanisation intensity.
4. Make a fine massage of your face, as shown in Fig. 3.

Recommendation

- Frequency of using the CLEANSING mode: 2-3 times a week (persons with more sensitive skin 1-2 times a week)
- Duration of treatment in the CLEANSING mode: 4 minutes (persons with more sensitive skin 1-2 minutes)

IMPORT

This mode is intended for easier penetration of beauty preparations into deeper layers of the skin. It combines negative ionisation, photon therapy using red light, heating and vibrations and allows efficient penetration of beauty preparations, such as moisturising and rejuvenating crèmes, hyaluronic acid, oligopeptides etc. through the skin.

Procedure (Fig. 4):

1. Apply the beauty preparation onto the skin.
2. Press **A1** to switch on the appliance, select the IMPORT mode and take the appliance in your hand, touching area **A6** with your hand/fingers.
3. Where necessary and considering sensitivity of the skin, use button **A2** to adjust temperature and **A3** to adjust galvanisation intensity.
4. Treat your face as shown in Fig. 4. As soon as you feel resistance on the application head, the preparation has been absorbed by the skin, and you can proceed on next part of your face.

Recommendation

- Frequency of using the IMPORT mode: 2-3 times a week (persons with more sensitive skin 1-2 times a week)
- Duration of treatment in the IMPORT mode: 4 minutes (persons with more sensitive skin 1-2 minutes)

LIFTING

This mode is intended for smoothening of face contours and making them finer. It combines, heating, ionisation, and EMS impulse therapy, providing effect on skin muscles and elastic collagen fibres.

Procedure (Fig. 5):

1. Apply the suitable jelly on application head **A4**.
2. Press **A1** to switch on the appliance, select the LIFTING mode and take the appliance in your hand, touching area **A6** with your hand/fingers.
3. Where necessary and considering sensitivity of your skin use button **A2** to adjust temperature and **A3** to adjust galvanisation intensity.
4. Treat your face as shown in Fig. 5, proceeding on individual parts of the face from the middle to the sides.

Recommendation

- Frequency of using the LIFTING mode: 1-2 times a week (persons with more sensitive skin once a week)
- Duration of treatment in the LIFTING mode: 4-8 minutes

EYE CARE

This mode is intended for skin care around the eyes. It combines heating, ionisation, vibrations and photon therapy using red light and allows easy penetration of the face mask and other crèmes through the skin.

Procedure (Fig. 6):

1. Apply the beauty preparation onto the skin, around your eyes.
2. Press **A1** to switch on the appliance, select the EYE CARE mode and take the appliance in your hand, touching area **A6** with your hand/fingers.
3. You can also switch on/off optional photon therapy using red light by brief pressing button **A3**, where necessary.
4. Make a fine massage of the face, around the eyes, as shown in Fig. 6.

Doporučení

- Frequency of using the EYE CARE mode: 2-3 times a week (persons with more sensitive skin 1-2 times a week)
- Duration of treatment in the EYE CARE mode: 4 minutes (persons with more sensitive skin 1-2 minutes)

COOL

This mode is intended for calming the skin. It combines cooling and photon therapy using blue light, which has anti-inflammatory effect, and reduces skin redness, acne problems and also allows easy penetration of calming preparations and rejuvenating face mask.

Procedure (Fig. 7):

1. Apply the required face mask.
2. Press **A1** to switch on the appliance and select the COOL mode.
3. Smoothly press on sensitive and/or red spots, as shown in Fig. 7.

Recommendation

- Frequency of using the COOL mode: 2-3 times a week
- Duration of treatment in the COOL mode: 3 minutes (maximum 5 seconds at each spot)

Automatic switching off

In CLEANSING, IMPORT, LIFTING and EYE CARE modes, the selected function will be deactivated after approx. 4 minutes (sound signal). After another approx. 1 minute, the appliance will be switched off (sound signal).

In COOL mode, the selected function will be deactivated after approx. 3 minutes (sound signal). After another approx. 1 minute, the appliance will be switched off (sound signal).

IV. MAINTENANCE

The appliance must be cleaned after every use! Use only soft cloth to wipe the surface of the appliance and the USB cable. Storage **D** can be rinsed with flowing water. Do not use any abrasive or aggressive cleaning agents!

Guarantee of 24 months does not apply to decrease in the accumulator's effective capacity caused by its use or age. Effective capacity gradually decreases depending on the method of use of the accumulator.

Expert service must do the large maintenance or to intervene into inside parts of the appliance!
Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!

V. ENVIRONMENT

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

Taking out the battery

Only take out the battery from the appliance after it is fully discharged, and after the end of the appliance's service life! In the course of operation, the battery does not need to be removed! For taking out the battery, remove rubber screw covers using a screwdriver, unscrew the screws and take out the battery.

VI. TECHNICAL DATA

Voltage (V)	Specified on the rating place
Appliance weight (g)	145
Appliance protection class	III.
Dimensions (H x D x W) in mm	51 x 163 x 58

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

The  symbol indicates a **WARNING**.

Do not use the appliance close to bath tubs, showers, bath basins or other water containers.

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© DATE 27/06/2022